English	2
Lithuanian	6
Danish	10
Estonian	14
Latvian	18
Norwegian	22
Finland	26
Swedish	30

Oclean Smart Sonic Electric
Toothbrush User Manual

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
Oclean Smart Sonic Electric Toothbrush User Manual

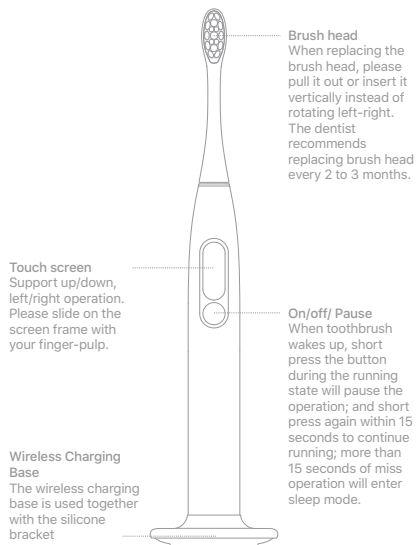Instructions

- First Time Use**
The toothbrush can be turned on by pressing the ON/OFF button for the first time on a brand-new toothbrush. The Bluetooth function will be activated and remain on. After installing the Oclean App, please complete the prompts to complete the connection and settings for the toothbrush.
- Install App & Toothbrush Settings**
You can use the app to customize your brushing programs - automating tooth cleaning, tooth whitening and gum massaging in one brushing session. There is no need to open the Oclean App each time you brush. Simply open the app every few days to sync the brushing data.
- Direct Use without App Settings**
You can also use the toothbrush without connecting to the app. It has a touchscreen display and an ON/OFF button. It offers touch control via the touch-sensitive screen.
- Brushing Plan and Performance Records**
Toothbrush will start default mode with "Standard Cleaning Plan" with 2 mins if no new customized brushing plan setting in APP.

After brushing, the display will show the effective brushing time, the brushing scores and the tips of the brushing area.

- How to Charge**
Please place the designated wireless charging base of the product on a smooth surface and connect the charging base to the plug. When power is applied, the indicator light on the base will flash. After mounting the silicone protective sleeve on the surface of the wireless charging base, place the toothbrush handle in the slot of the wireless charging base. If not charged, the indicator light will not lit. When there is a foreign body on the wireless charger or when it is not being charged properly, the indicator light will flash quickly. When the toothbrush is being charged properly, the indicator light will turn on a breathing-like display effect and the battery will generally be fully charged in around 3.5 hours.

Follow the Oclean Official Facebook Page

Oclean Facebook Page
www.facebook.com/ocleanglobal

Search "Oclean" on App Store or Google Play to download

Product specifications

Product Name	Oclean Smart Sonic Electric Toothbrush
Model	Y2087
Version	Pro Elite
Emission power	0 dBm
IP Rating	IPX7
Power of the toothbrush body	1.5W
Charging	Oclean USB charging base
Input	5V = 0.3A
Operating time*	20 days (level 32 of intensity) - 40 days (level 1 of intensity)
Battery Capacity	800 mAh
Charging time	3.5h
Implementation standard	GB 4706.1-2005 GB 4706.59-2008
Bluetooth	Bluetooth 5.0 above

* Under the cleaning mode, it is calculated at 2 minutes/time, 2 times/day working condition.

R-R-y1T-Y2087

IPX7

210-155489

Remarks

- This indicates that the content of toxic and hazardous substances in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011.
- × This indicates that the content of toxic and hazardous substances in at least one of the homogeneous materials of this component is more than the limited quantity as specified in GB/T 26572-2011.
- * Except for very few materials, due to technical development reasons, most of the materials used are environmentally friendly. Please use this product with confidence during the environmental protection period.

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Warranty Card

- These clauses only apply to Oclean products sold in the country or region where the purchase is made, if there is any inconsistency or omission with the national policy, the national policy shall prevail.
- For performance issues caused by the product itself, one may avail service as per below regulations with valid proof of purchase:
 - A customer can avail warranty service for free within 24 months from the date of purchase if there is a performance failure caused by the product itself.
 - In order to avail replacement and warranty services, one should show the valid proof of purchase.
 - Other product related materials (such as Packaging box, Product Manual etc.) are not covered in warranty.

Product model:	Serial Number:					
	Name:					
User Information	Tel:					
	Email:					
	Address:					
	Date of purchase:					
	Invoice number:					
Problem Description:						

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Follow the Oclican Official
Facebook Page

Oclican Facebook Page
www.facebook.com/oclicanglobal

Search "Oclican" on App Store or Google Play to download

- Product specifications

Product Name	Oclican Smart Sonic Electric Toothbrush
Model	Y2087
Version	Pro Elite
Emission power	0 dBm
IP Rating	IPX7
Power of the toothbrush body	1.5W
Charging	Oclican USB charging base
Input	5V = 0.3A
Operating time*	20 days (level 32 of intensity) - 40 days (level 1 of intensity)
Battery Capacity	800 mAh
Charging time	3.5h
Implementation standard	GB 4706.1-2005 GB 4706.59-2008
Bluetooth	Bluetooth 5.0 above

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210-155489

Remarks

- This indicates that the content of toxic and hazardous substances in all homogeneous materials of this component is less than the limited quantity as specified in GB/T 26572-2011.
- ✗ This indicates that the content of toxic and hazardous substances in at least one of the homogeneous materials of this component is more than the limited quantity as specified in GB/T 26572-2011.
- * Except for very few materials, due to technical development reasons, most of the materials used are environmentally friendly. Please use this product with confidence during the environmental protection period.

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